



TD52 Load Cell Cable Adaptor Kit


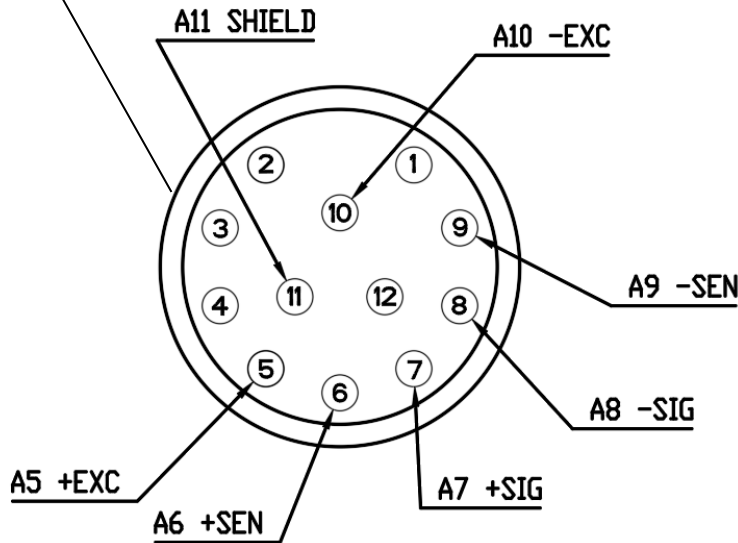


Table1

Pin	Connection
A5	+EXC
A6	+SEN
A7	+SIG
A8	-SIG
A9	-SEN
A10	-EXC
A11	SHIELD

Figure 1

- 1** Open the Load Cell Cable adaptor. The adaptor contains head, body and cap.
- 2** Pass the Load Cell Cable through the adaptor's cap and body.
- 3** Weld the Load Cell Cable carefully to the adaptor's head one by one according to table 1 and figure 1.

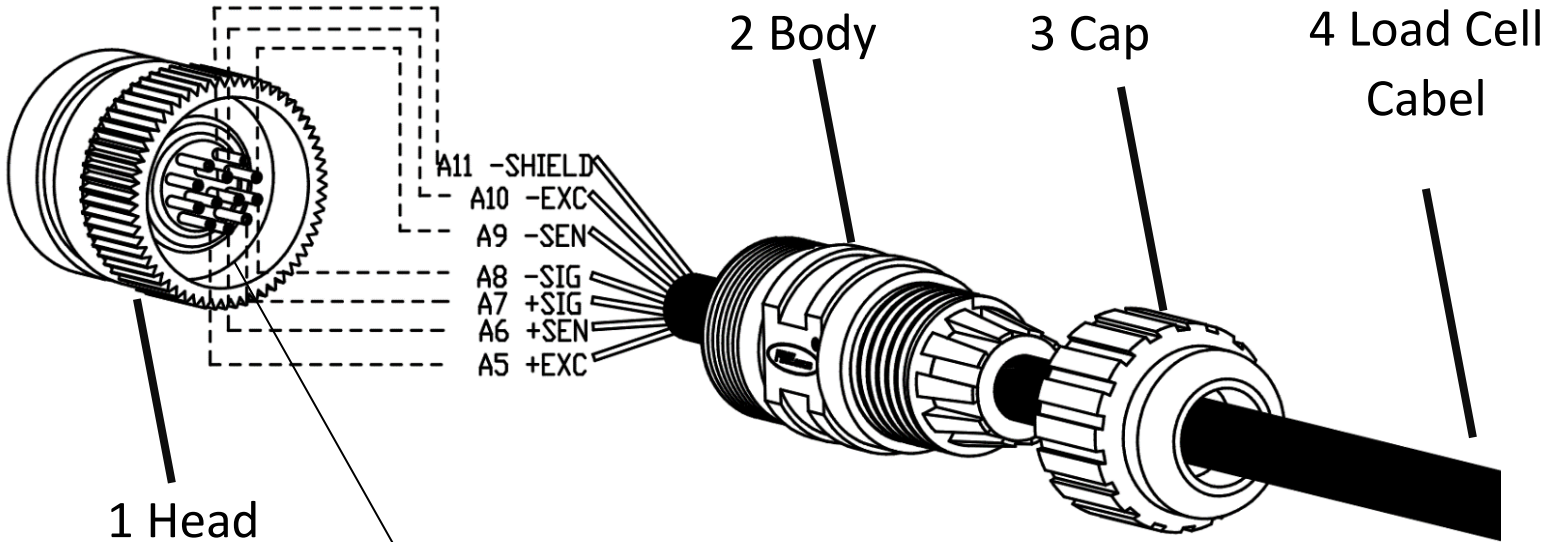

4 Close the Load Cell Cable adaptor. The following figure 2 shows the final status.

Note: the Load Cell Cable adaptor is not waterproof.

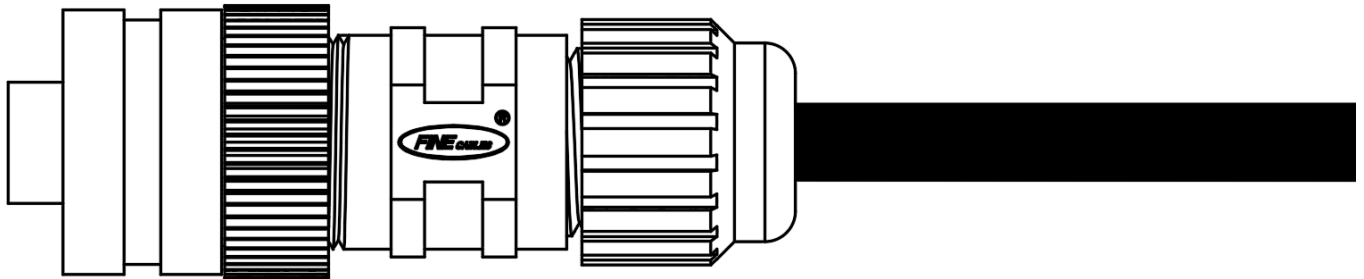


Figure 2



P/N 30553821 A ©Ohaus Corporation 2019, all right reserved.

Printed in China